

DPEX X

Precision Beyond Measure



APEX LOCATOR

Comprehensive Configuration

Not autoclavable

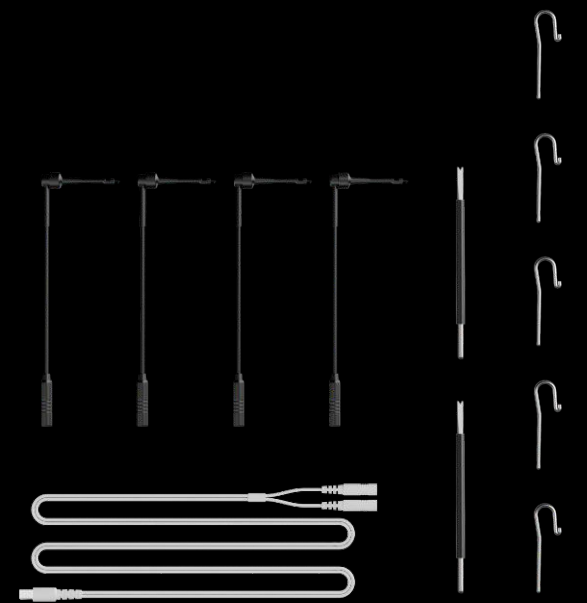
Measuring wire x 1 (Soft texture)

Autoclavable

File clip x 4 (Coated with biological silica gel)

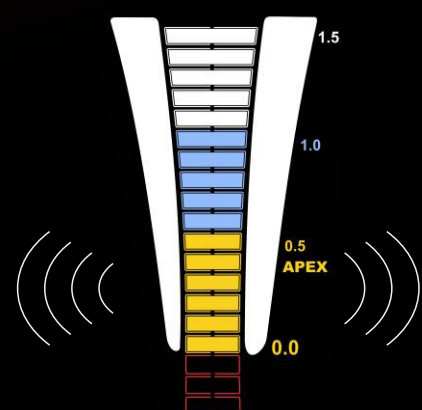
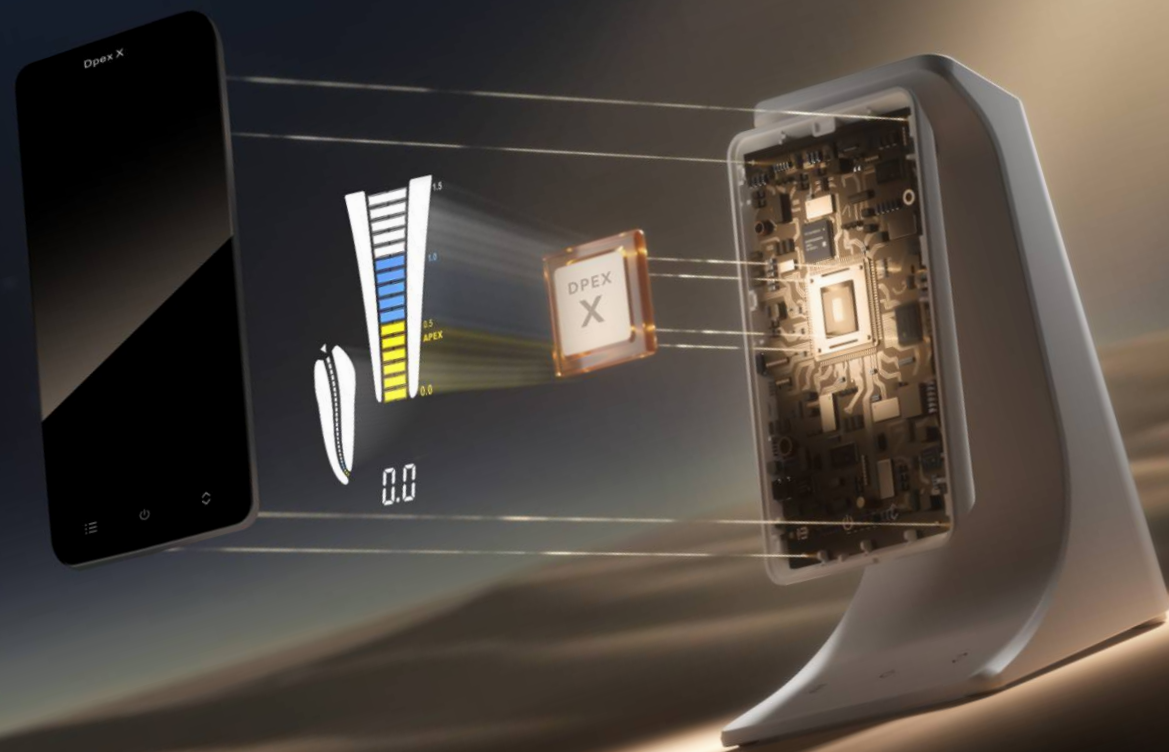
Lip hook x 5 (Made of medical-grade stainless steel)

Apex locator probe x 2



Guilin Woodpecker Medical Instrument Co., Ltd.
 Information Industrial Park, Guilin National High-Tech
 Zone, Guilin, Guangxi, 541004 P. R. China
 Sales Dept.: +86-0773-5873196/2350599
 After-sales Service Dept.: +86-0773-5827898
 E-mail: woodpecker4@glwoodpecker.com
 Website: http://www.glwoodpecker.com
 ZMN-XC-169 V1.0

DISTRIBUTED BY:

New Circuit Design
Minimized InterferenceSound Recognition
Capable Aid in Operation

Two sound modes
 Adjust volume and mode as you wish



Working Time: 7h Charging Time: 3h Battery: 3.7V/2000mAh Display: 3.8-inch LCD screen

Large Color Screen
Clearer and Visual Display

3.8-inch LCD screen

Magnification to display the apical third

Clinical recommendation:

Measure to the 0.0 position and manually
 subtract 0.5~1mm for the final working length

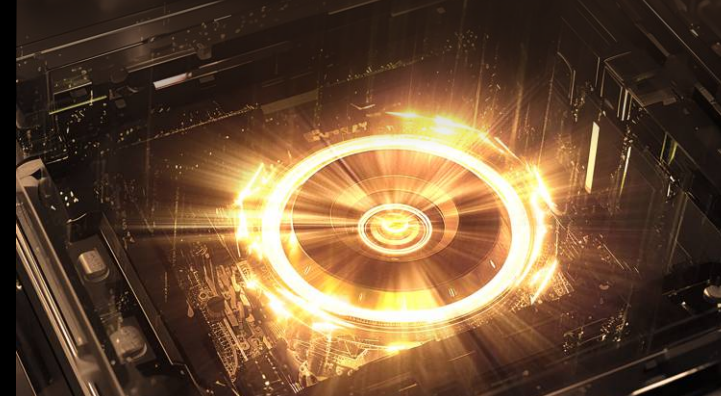


Four Working Modes

6.2s to reach the maximum pulse
 intensity at the quickest

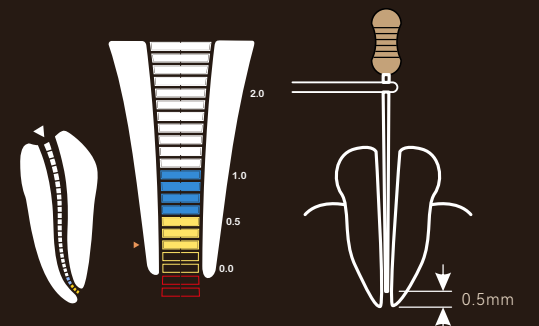
Mode 1: 42.5s Mode 2: 32s

Mode 3: 24s Mode 4: 6.2s



Support Apical Stop Reset

The apical stop can be freely set according
 to clinical needs.
 Reset the display point that is 0.5mm from the apex.
 Provide warning prompts based on the
 preset apical stop.

Added Sleep Mode
Automatic Wake-up

Three sleep durations available

The device will be woken up upon the
 measurement in sleep mode

10s 20s 30s

